**NZYtech** 

CZ0364\_UG\_EN\_V2302

# Carbohydrate Binding Module 13A, Ruminococcus flavefaciens

# (CBM13)


Catalogue number CZ03641 CZ03642 Presentation 1 mg 3 x 1 mg

## Description

Carbohydrate Binding Module 13A (CBM13) is a Carbohydrate Binding Protein originating from *Ruminococcus flavefaciens*. The recombinant CBM13, purified from *Escherichia coli*, is a single-domain protein belonging to the Carbohydrate Binding Module family 13 (CBM13, see more details at <u>www.cazy.org</u>). The protein is supplied in a solution containing 35 mM NaHepes buffer (pH 7.5), 750 mM NaCl, 200 mM Imidazole, 3.5 mM CaCl<sub>2</sub>, and 25% (v/v) glycerol, at a concentration of 1 mg/mL. Bulk quantities of this product can be made available upon request. To place an order, simply visit our website. We offer fast and secure shipping worldwide.

## **Electrophoretic Purity**

The molecular integrity and purity of CBM13 were evaluated using sodium dodecyl sulfate-polyacrylamide gel electrophoresis (SDS-PAGE), followed by BlueSafe staining (MB15201) (Figure 1).

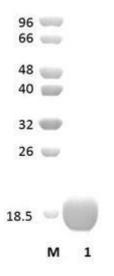

**Figure 1**. SDS-PAGE analysis of CBM13 was conducted in (Lane 1), employing a 14% polyacrylamide gel. The enzyme exhibits a band corresponding to a molecular weight of approximately 17,05 kDa. Lane M contains a Protein Marker for reference.

#### Storage temperature

The protein should be stored at -30°C to -15°C in a constant temperature freezer. The protein will remain stable till the expiry date if stored as specified.

#### **Ligand specificity**

CBM13 binds to arabinoxylan.

# **Customer Support**

Our dedicated customer support team is always ready to assist you with any questions or technical issues you may have. Reach us via email at info@nzytech.com.

# **Quality control assay**

Protein purity is determined to be ≥50%, as assessed by SDS-PAGE and subsequent BlueSafe staining (MB15201).

For life science research only. Not for use in diagnostic procedures.

NZYtech Lda. Estrada do Paço do Lumiar, Campus do Lumiar - Edifício E, R/C, 1649-038 Lisboa, Portugal Tel.:+351.213643514 Fax: +351.217151168 www.nzytech.com