B NZYtech

CZ0469_UG_EN_V2302

Carbohydrate Binding Module 79A, Ruminococcus flavefaciens

(CBM79)

 Catalogue number
 Pro

 CZ04691
 0.5

 CZ04692
 3 >

Presentation 0.5 mg 3 x 0.5 mg

Description

Carbohydrate Binding Module 79A (CBM79) is a Carbohydrate Binding Protein originating from *Ruminococcus flavefaciens*. The recombinant CBM79, purified from *Escherichia coli*, is a single-domain protein belonging to the Carbohydrate Binding Module family 79 (CBM79, see more details at <u>www.cazy.org</u>). The protein is supplied in a solution containing 35 mM NaHepes buffer (pH 7.5), 750 mM NaCl, 200 mM Imidazole, 3.5 mM CaCl₂, and 25% (v/v) glycerol, at a concentration of 0.5 mg/mL. Bulk quantities of this product can be made available upon request. To place an order, simply visit our website. We offer fast and secure shipping worldwide.

Electrophoretic Purity

The molecular integrity and purity of CBM79 were evaluated using sodium dodecyl sulfate-polyacrylamide gel electrophoresis (SDS-PAGE), followed by BlueSafe staining (MB15201) (Figure 1).

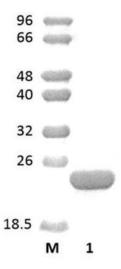


Figure 1. SDS-PAGE analysis of CBM79 was conducted in (Lane 1), employing a 14% polyacrylamide gel. The enzyme exhibits a band corresponding to a molecular weight of approximately 20,97 kDa. Lane M contains a Protein Marker for reference.

Storage temperature

The protein should be stored at -30°C to -15°C in a constant temperature freezer. The protein will remain stable till the expiry date if stored as specified.

Ligand specificity

CBM79 binds to 1,3-1,4-β-glucans, soluble 1,4-β-glucans and xyloglucan. The biochemical properties of CBM79 are detailed in the referenced publication(s) provided below.

Reference

Venditto et al. (2016) Proc Natl Acad Sci U S A. 113(26):7136-41.

PDB/3D code: 4D3L, 4V1K, 4V1L.

Customer Support

Our dedicated customer support team is always ready to assist you with any questions or technical issues you may have. Reach us via email at info@nzytech.com.

Quality control assay

Protein purity is determined to be ≥90%, as assessed by SDS-PAGE and subsequent BlueSafe staining (MB15201).

For life science research only. Not for use in diagnostic procedures.

NZYtech Lda. Estrada do Paço do Lumiar, Campus do Lumiar - Edifício E, R/C, 1649-038 Lisboa, Portugal Tel.:+351.213643514 Fax: +351.217151168 www.nzytech.com