

CZ0803_UG_EN_V2302

Aceric acid hydrolase 127A, Bacteroides thetaiotaomicron

*Bt*Ach127A (GH127)

Catalogue number	Presentation
CZ08031	1 mg
CZ08032	3 x 1 mg

Description

Aceric acid hydrolase 127A (*Bt*Ach127A), assigned the E.C. number 3.2.1.-, is a derivative of *Bacteroides thetaiotaomicron*. It is a 3-C-carboxy-5-deoxy-L-xylose (aceric acid) hydrolase. The recombinant *Bt*Ach127A, purified from *Escherichia coli*, is a single-domain Glycoside Hydrolase family 127 (GH127) enzyme (see more details at <u>www.cazy.org</u>). The protein is supplied in a solution containing 35 mM NaHepes buffer (pH 7.5), 750 mM NaCl, 200 mM Imidazole, 3.5 mM CaCl₂, and 25% (v/v) glycerol, at a concentration of 1 mg/mL. Bulk quantities of this product can be made available upon request. To place an order, simply visit our website. We offer fast and secure shipping worldwide.

Electrophoretic Purity

The molecular integrity and purity of *Bt*Ach127A (GH127) were evaluated using sodium dodecyl sulfate-polyacrylamide gel electrophoresis (SDS-PAGE), followed by BlueSafe staining (MB15201) (Figure 1).

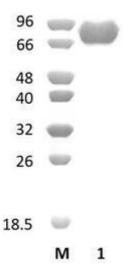


Figure 1. SDS-PAGE analysis of *Bt*Ach127A (GH127) was conducted in (Lane 1), employing a 14% polyacrylamide gel. The enzyme exhibits a band corresponding to a molecular weight of approximately 79,80 kDa. Lane M contains a Protein Marker for reference.

Storage temperature

The protein should be stored at -30°C to -15°C in a constant temperature freezer. The protein will remain stable till the expiry date if stored as specified.

Substrate specificity

BtAch127A (GH127) hydrolyses the L-Acef- α 1,3-L-Rhap linkage in rhamnogalacturonan II (RGII).

Temperature and pH optima

The enzyme exhibits optimal activity within a pH range of 6.5-7.5 and at a temperature of 37°C. Maximal enzymatic activity is achieved at pH 7 and a consistent temperature of 37°C.

Enzyme activity

The substrate specificity and kinetic properties of *Bt*Ach127A (GH127) are detailed in the referenced publication provided below. To perform enzyme assays and determine specific activity values, adhere to the methodology outlined in the cited paper(s).

Reference

Ndeh et al. (2017) Nature. 544(7648):65-70.

Mahowald et al. (2009) Proc Natl Acad Sci U S A. 106(14):5859-5864.

Customer Support

Our dedicated customer support team is always ready to assist you with any questions or technical issues you may have. Reach us via email at info@nzytech.com.

Quality control assay

Protein purity is determined to be ≥90%, as assessed by SDS-PAGE and subsequent BlueSafe staining (MB15201).

For life science research only. Not for use in diagnostic procedures.

NZYtech Lda. Estrada do Paço do Lumiar, Campus do Lumiar - Edifício E, R/C, 1649-038 Lisboa, Portugal Tel.:+351.213643514 Fax: +351.217151168 www.nzytech.com