

CZ0890 UG EN V2302

Carbohydrate Binding Module 16A, Caldanaerobius polysaccharolyticus (CBM16)

Catalogue number Presentation

CZ08901 1 mg CZ08902 3 x 1 mg

Description

Carbohydrate Binding Module 16A (CBM16) is a Carbohydrate Binding Protein originating from *Caldanaerobius polysaccharolyticus*. The recombinant CBM16, purified from *Escherichia coli*, is a single-domain protein belonging to the Carbohydrate Binding Module family 16 (CBM16, see more details at www.cazy.org). The protein is supplied in a solution containing 35 mM NaHepes buffer (pH 7.5), 750 mM NaCl, 200 mM Imidazole, 3.5 mM CaCl₂, and 25% (v/v) glycerol, at a concentration of 1 mg/mL. Bulk quantities of this product can be made available upon request. To place an order, simply visit our website. We offer fast and secure shipping worldwide.

Electrophoretic Purity

The molecular integrity and purity of CBM16 were evaluated using sodium dodecyl sulfate-polyacrylamide gel electrophoresis (SDS-PAGE), followed by BlueSafe staining (MB15201) (Figure 1).

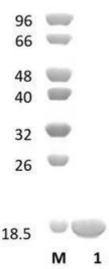


Figure 1. SDS-PAGE analysis of CBM16 was conducted in (Lane 1), employing a 14% polyacrylamide gel. The enzyme exhibits a band corresponding to a molecular weight of approximately 18,17 kDa. Lane M contains a Protein Marker for reference.

Storage temperature

The protein should be stored at -30°C to -15°C in a constant temperature freezer. The protein will remain stable till the expiry date if stored as specified.

Ligand specificity

CBM16 binds to glucose-, mannose-, and glucose/mannose-configured substrates. The biochemical properties of CBM16 are detailed in the referenced publication(s) provided below.

Reference

Bae et al. (2008) J Biol Chem. 283(18):12415-25.

Su et al. (2010) J Biol Chem. 285(45):34665-76.

PDB/3D code: 2ZEW, 2ZEX, 2ZEY, 2ZEZ.

Customer Support

Our dedicated customer support team is always ready to assist you with any questions or technical issues you may have. Reach us via email at info@nzytech.com.

Quality control assay

Protein purity is determined to be ≥90%, as assessed by SDS-PAGE and subsequent BlueSafe staining (MB15201).

For life science research only. Not for use in diagnostic procedures.