

CZ1038\_UG\_EN\_V2302

# Carbohydrate Binding Module 20A, Bacillus sp.

# (GFP-CBM20)


Catalogue number	Presentation	
CZ10381	1 mg	
CZ10382	3 x 1 mg	

# Description

Carbohydrate Binding Module 20A (GFP-CBM20) is a Carbohydrate Binding Protein originating from *Bacillus sp.*. The recombinant GFP-CBM20, purified from *Escherichia coli*, is a modular protein belonging to the Carbohydrate Binding Module family 20 (CBM20, see more details at <u>www.cazy.org</u>) fused to an N-terminal Green Fluorescent Protein (GFP). This GFP protein derivative is particularly recommended for subcellular localization studies, which allows for real-time tracking and imaging in living cells. The protein is supplied in a solution containing 35 mM NaHepes buffer (pH 7.5), 750 mM NaCl, 200 mM Imidazole, 3.5 mM CaCl<sub>2</sub> and 3.2 M ammonium sulphate, at a concentration of 1 mg/mL. Bulk quantities of this product can be made available upon request. To place an order, simply visit our website. We offer fast and secure shipping worldwide.

# **Electrophoretic Purity**

The molecular integrity and purity of GFP-CBM20 were evaluated using sodium dodecyl sulfate-polyacrylamide gel electrophoresis (SDS-PAGE), followed by BlueSafe staining (MB15201) (Figure 1).

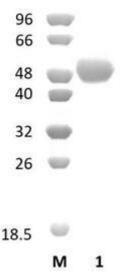

**Figure 1**. SDS-PAGE analysis of GFP-CBM20 was conducted in (Lane 1), employing a 14% polyacrylamide gel. The enzyme exhibits a band corresponding to a molecular weight of approximately 49,25 kDa. Lane M contains a Protein Marker for reference.

#### Storage temperature

The protein should be stored at 2°C to 8°C in a constant temperature freezer. The protein will remain stable till the expiry date if stored as specified.

#### **Ligand specificity**

GFP-CBM20 binds to cyclodextrins. The biochemical properties of GFP-CBM20 are detailed in the referenced publication(s) provided below.

# **Assay conditions**

For optimal recovery of GFP-CBM20 activity, carry out the following procedure: centrifuge the necessary volume of the precipitated protein suspension at 13,000 x g for a duration of 5 minutes. Subsequently, decant the ammonium sulphate supernatant and resuspend the resultant pellet in an equivalent volume of solution, comprising 20 mM Tris-HCl (pH 7.5), 20 mM NaCl, and 5 mM CaCl2. Following resuspension, proceed to the appropriate assay as dictated by your experimental requirements.

# Reference

Kimura et al. (1990) Agric Biol Chem. 54(3):641-8.

# **Customer Support**

Our dedicated customer support team is always ready to assist you with any questions or technical issues you may have. Reach us via email at info@nzytech.com.

# **Quality control assay**

Protein purity is determined to be ≥75%, as assessed by SDS-PAGE and subsequent BlueSafe staining (MB15201).

For life science research only. Not for use in diagnostic procedures.

NZYtech Lda. Estrada do Paço do Lumiar, Campus do Lumiar - Edifício E, R/C, 1649-038 Lisboa, Portugal Tel.:+351.213643514 Fax: +351.217151168 www.nzytech.com