

TX1B3 PDZ Domain, *Homo sapiens*

Catalogue number	Presentation
PD01201	0.25 mg
PD01202	3 x 0.25 mg

Description

TX1B3 PDZ Domain from *Homo sapiens* is a recombinant protein purified from *Escherichia coli*. The protein is provided in 35 mM NaHepes buffer, pH 7.5, 750 mM NaCl, 200 mM imidazol, 3.5 mM CaCl₂ and 25% (v/v) glycerol, at a 0.5 mg/mL concentration. Bulk quantities of this product can be made available upon request. To place an order, simply visit our website. We offer fast and secure shipping worldwide.

Electrophoretic Purity

TX1B3 PDZ Domain purity was evaluated using sodium dodecyl sulfate-polyacrylamide gel electrophoresis (SDS-PAGE), followed by BlueSafe staining (MB15201) (Figure 1).

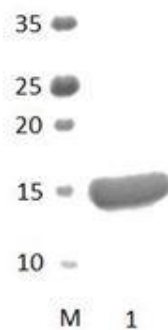


Figure 1. SDS-PAGE analysis of TX1B3 PDZ Domain was conducted in (Lane 1), employing a 14% polyacrylamide gel. The enzyme exhibits a band corresponding to a molecular weight of approximately 16.03 kDa. Lane M contains a Protein Marker for reference.

Storage temperature

The protein should be stored at -30°C to -15°C in a constant temperature freezer. The protein will remain stable till the expiry date if stored as specified.

Protein Sequence

QPVTAVVQRVEIHKLRQGENLILGFSIGGGIDQDPSQNPFSQDKTKGIYVTRVSEGGPAEIAQLQIGDKIMQVNGWDMTMMVTHDQARKRLTKRSEEVVRLLVTRQS
LQKAVQQSMLS

Number of PDZ in native protein

1

Temperature and pH optima

The protein exhibits optimal activity within a pH of 7 and at a temperature of 36.5 °C.

PDB code

2KG2, 4NNM, 3SFJ, 4NNL, 4E3B

Reference

J. Biomol. NMR 45 (3), 329-334 (2009); J. Biol. Chem. 288 (7), 5114-5126 (2013)

Customer Support

Our dedicated customer support team is always ready to assist you with any questions or technical issues you may have. Reach us via email at info@nzytech.com.

Quality control assay

Protein purity is determined to be $\geq 90\%$, as assessed by SDS-PAGE and subsequent BlueSafe staining (MB15201).

For life science research only. Not for use in diagnostic procedures.

NZYtech Lda. Estrada do Paço do Lumiar, Campus do Lumiar - Edifício E, R/C, 1649-038 Lisboa, Portugal Tel.:+351.213643514 Fax:
+351.217151168 www.nzytech.com